

Reliable, Affordable 58 Litre Toploading Autoclave!

Only \$8,750

Was \$14,726

Code: TOMYSX-500E

includes 2 stainless steel mesh baskets




Five sterilization profiles to select from (in addition to normal sterilization)

Three operating temperature ranges:

Warming 45 to 95°C, Heating 45 to 104°C and Sterilizing 105 to 135°C

- Integrated (removable) condensate bottle to trap condensate from chamber during venting
- Requires no plumbing so can be moved to wherever it's needed
- Powered from a 240v 50hz 15amp supply
- Simple lid opening using either hand or foot
- Basket Dimensions: 300mm Ø x 181mm H
- Grade 304 stainless steel for easy cleaning and durability
- Variable speed exhaust
- Safety devices include: water level sensor, current leakage breaker, lid interlock, overheat prevention, over pressure prevention, open temperature sensor detection & safety valve
- Digital temperature & time display with analogue gauge for pressure
- 2-year warranty

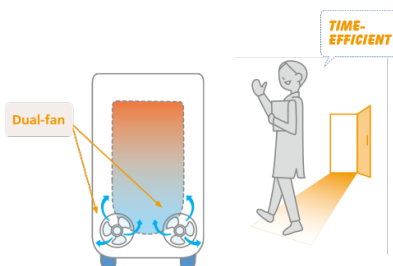

The SX-E series Autoclaves are standard-equipped with a dual-fan cooling system, shortening overall operation time.

Good range of accessories available including:

- Punched metal stainless steel baskets
- Solid stainless steel baskets
- Hybrid basket: solid lower third with punched metal stainless steel upper
- Durham tube stand
- Printer
- External sensor (wandering probe)
- Preventative maintenance contracts

Prices and discounts may be subject to change and only available while stocks last. Not valid with any other special offer or discounts and cannot be applied in retrospect. All prices in local currency excluding GST and local delivery charges. Products advertised may be for research use only (RUO) or in vitro diagnostic use (IVD). Please refer to manufacturer's product insert for further information. For product use outside of that specified on the manufacturer's labelling and/or instructions for use, Bio-Strategy does not take responsibility.

Ref. Code:
TOMY-SX-500E

VALID UNTIL:
30 Apr 2024