

## Separator Kit For The Disposal Of Petri Dishes

The disposable petri dish contains both bio-waste (plastic) and media with cultured bacteria that is potentially infectious.


The Separator Kit not only separates bio-waste from media but is the easiest solution for sterilizing the bio-waste and media in a single autoclave cycle.

This accessory works for disposable (plastic) petri dishes only. Not suitable for glass petri dishes.



Melted petri dishes remain in the long basket after the sterilisation.

The agar melts and falls down into the bucket.

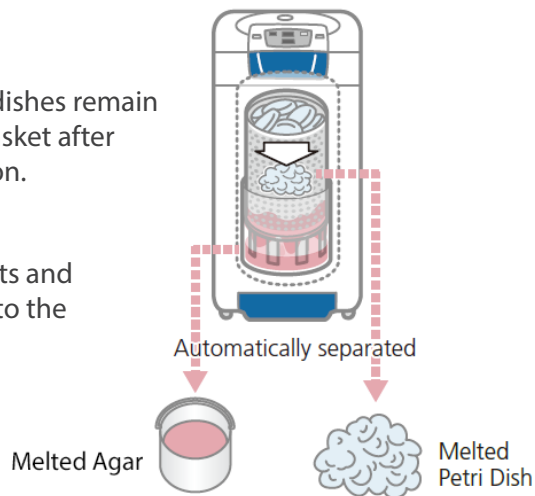

Watch >




## ST-ZERO Bio waste Sterilisation Container

Bio waste disposables (infectious sample) should be sterilized before removing from the user's facility.


How do users carry these samples from bio safety cabinet to Autoclave? ST-ZERO is one solution for this issue and can be used with all models of TOMY Autoclaves.

Designed for liquid containment with flat placement and vertical mounting. Air flow will not be disturbed inside the work area of bio-safety cabinet.

- easy installation
- capacity for up to 25 x 25ml disposable pipettes
- leak proof while sterilising in vertical autoclaves

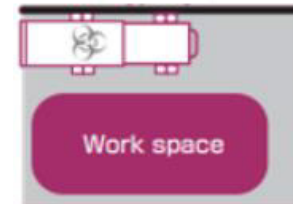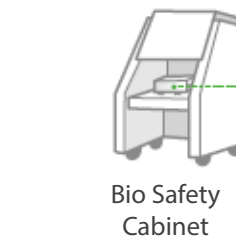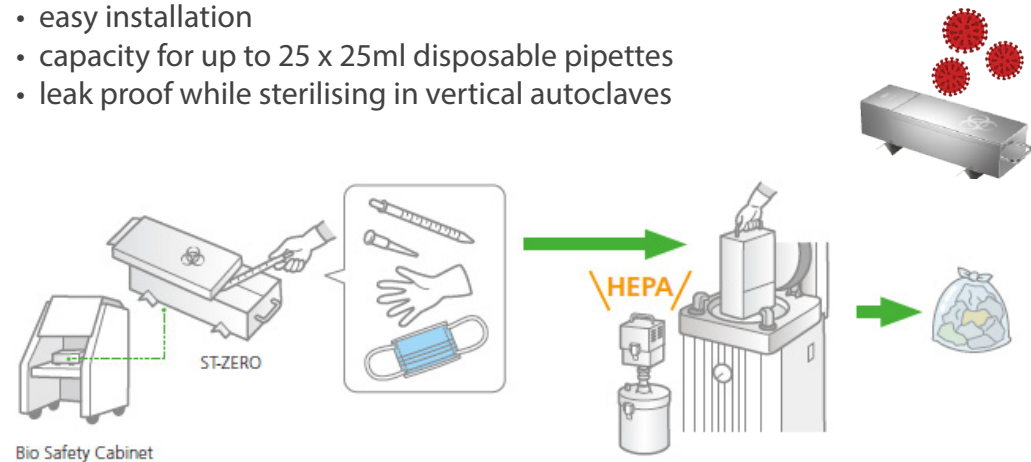

Installed horizontally

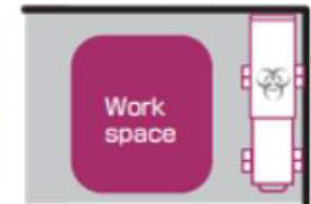

Installed vertically