

## Askion ColdEye+ Scanner

### The only scanner for all applications.


#### Performance features:

- Sample recognition within less than 5 seconds
- Simultaneous identification of several samples in SBS rack (96, 48, 24, etc.) or box format (9x9, 10x10, etc.)
- Reliable sample recognition at ambient temperatures from +25°C to -150°C
- Reads vial bottom code and rack / box barcode (1D/2D)
- Available as a portable standalone version or as a retrofittable installation solution for the workbench
- Sample format independent
- Integrated analysis software (for portable version)
- Integrated power supply for up to 5 h (at -150°C, for portable version)
- Enables best possible sample quality, due to minimal temperature changes during the scanning process

The ColdEye+ is the first and only portable scanner worldwide that enables safe and reliable sample identification at warm and cryogenic temperatures. This scanner works without impairment at room temperature (+ 25°C), in dry ice (approx. -80°C) up to the vapor phase of liquid nitrogen (up to -150°C). In addition, the scanner is independent of the sample format, which means that it can be used for a wide variety of laboratory applications.

By using ColdEye+, the process in the laboratory (e.g. inventory of -20°C or -80°C stores) can become safer (permanently identical temperatures for the samples) and more efficient (fast scan sequence of several samples / SBS racks/ Boxes one after the other). ColdEye+ is a modern scanner that simplifies everyday laboratory work for the user while maintaining sample quality.

**Contact  
Bio-Strategy for  
more details!**

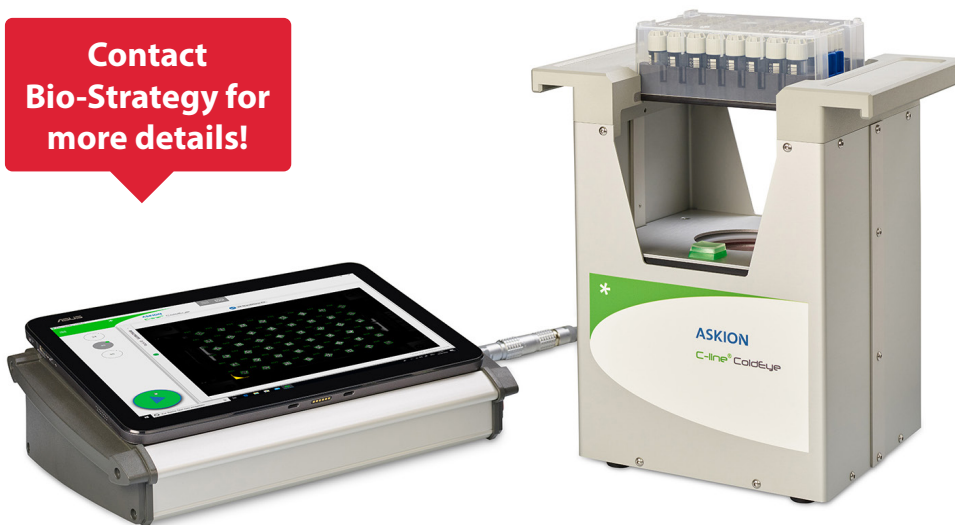

Contact us for more information  
and to place your orders

**Bio-Strategy Pty Ltd.**

T: 1800 008 453

E: [sales.au@bio-strategy.com](mailto:sales.au@bio-strategy.com)

[www.bio-strategy.com](http://www.bio-strategy.com)

[shop.bio-strategy.com](http://shop.bio-strategy.com)

 [linkedin.com/company/bio-strategy](https://www.linkedin.com/company/bio-strategy)

Rev. April 2021

Prices and discounts may be subject to change and only available while stocks last. Not valid with any other special offer or discounts and cannot be applied in retrospect. All prices in AUD excluding GST and local delivery charges. Products advertised may be for research use only (RUO) or in vitro diagnostic use (IVD). Please refer to manufacturer's product insert for further information. For product use outside of that specified on the manufacturer's labelling and/or instructions for use, Bio-Strategy does not take responsibility.